

Chromogranin A Polyclonal Antibody

Catalog Number:E-AB-40339



Note: Centrifuge before opening to ensure complete recovery of vial contents.


Description

Reactivity	Human,Mouse
Immunogen	Recombinant Mouse Chromogranin-A protein
Host	Rabbit
Isotype	IgG
Purification	Affinity purification
Conjugation	Unconjugated
Formulation	PBS with 0.05% Proclin300, 50% glycerol, pH7.3.

Applications Recommended Dilution


WB	1:500-1:1000
IHC	1:100-1:300

Data

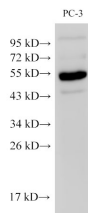


Western Blot analysis of PC-3 cells using Chromogranin A Polyclonal Antibody at dilution of 1:500


Observed Mw:51 kDa
Calculated Mw:51 kDa

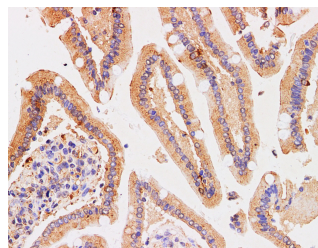

Immunohistochemistry of paraffin-embedded Human duodenum using Chromogranin A Polyclonal Antibody at dilution of 1:200

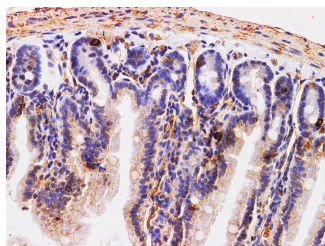

Immunohistochemistry of paraffin-embedded Mouse small intestine using Chromogranin A Polyclonal Antibody at dilution of 1:200

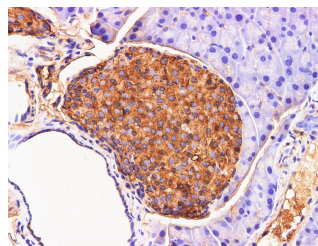

Immunohistochemistry of paraffin-embedded Mouse pancreas using Chromogranin A Polyclonal Antibody at dilution of 1:200

Preparation & Storage

Storage Store at -20°C. Avoid freeze / thaw cycles.

Background

Chromogranin A is a member of the granin family of neuroendocrine secretory proteins. It is located in secretory vesicles of neurons and endocrine cells. Chromogranin A is the precursor to several functional peptides including

For Research Use Only

A Reliable Research Partner in Life Science and Medicine

Tel: 400-999-2100

Email: techsupport@elabscience.cn

Web: www.elabscience.cn

Chromogranin A Polyclonal Antibody

Catalog Number:E-AB-40339



vasostatin,pancreastatin,catestatin and parastatin. These peptides negatively modulate the neuroendocrine function of the releasing cell (autocrine) or nearby cells (paracrine). CgA is one of the most used tumor markers in NET's (neuroendocrine tumors) ,and elevated CgA concentrations have been demonstrated in serum or plasma of patients with different types of these tumors.

For Research Use Only

A Reliable Research Partner in Life Science and Medicine

Tel: 400-999-2100

Email: techsupport@elabscience.cn

Web: www.elabscience.cn