

ELP1 Polyclonal Antibody

Catalog Number:E-AB-52014



Note: Centrifuge before opening to ensure complete recovery of vial contents.


Description

Reactivity	Human
Immunogen	Synthetic peptide of human ELP1
Host	Rabbit
Isotype	IgG
Purification	Affinity purification
Conjugation	Unconjugated
Formulation	PBS with 0.05% NaN ₃ and 40% Glycerol,pH7.4


Applications Recommended Dilution

IHC	1:30-1:150
ELISA	1:5000-1:10000

Data

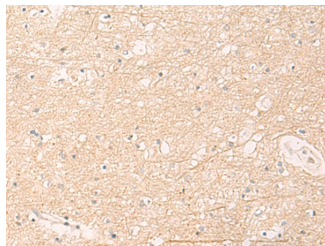

Immunohistochemistry of paraffin-embedded Human brain tissue using ELP1 Polyclonal Antibody at dilution of 1:30(×200)

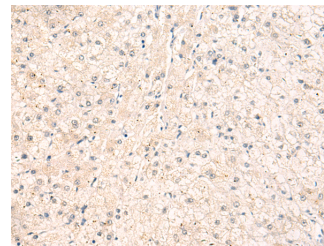

Immunohistochemistry of paraffin-embedded Human liver cancer tissue using ELP1 Polyclonal Antibody at dilution of 1:30(×200)

Preparation & Storage

Storage Store at -20°C. Avoid freeze / thaw cycles.

Background

The protein encoded by this gene is a scaffold protein and a regulator for three different kinases involved in proinflammatory signaling. The encoded protein can bind NF-kappa-B-inducing kinase and I-kappa-B kinases through separate domains and assemble them into an active kinase complex. Mutations in this gene have been associated with familial dysautonomia. Alternative splicing results in multiple transcript variants encoding different isoforms.

For Research Use Only

A Reliable Research Partner in Life Science and Medicine

Tel: 400-999-2100

Email: techsupport@elabscience.cn

Web: www.elabscience.cn