

CDKN1B Polyclonal Antibody

Catalog Number:E-AB-52142



Note: Centrifuge before opening to ensure complete recovery of vial contents.


Description

Reactivity	Human, Mouse
Immunogen	Fusion protein of human CDKN1B
Host	Rabbit
Isotype	IgG
Purification	Affinity purification
Conjugation	Unconjugated
Formulation	PBS with 0.05% NaN ₃ and 40% Glycerol,pH7.4


Applications Recommended Dilution

IHC	1:50-1:300
ELISA	1:5000-1:10000

Data

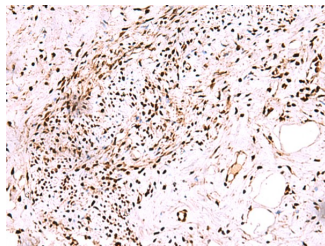

Immunohistochemistry of paraffin-embedded Human cervical cancer tissue using CDKN1B Polyclonal Antibody at dilution of 1:75(×200)

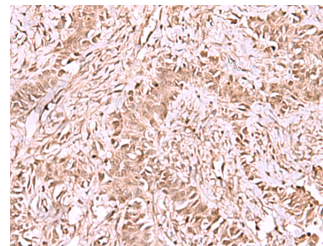

Immunohistochemistry of paraffin-embedded Human ovarian cancer tissue using CDKN1B Polyclonal Antibody at dilution of 1:75(×200)

Preparation & Storage

Storage Store at -20°C. Avoid freeze / thaw cycles.

Background

This gene encodes a cyclin-dependent kinase inhibitor, which shares a limited similarity with CDK inhibitor CDKN1A/p21. The encoded protein binds to and prevents the activation of cyclin E-CDK2 or cyclin D-CDK4 complexes, and thus controls the cell cycle progression at G1. The degradation of this protein, which is triggered by its CDK dependent phosphorylation and subsequent ubiquitination by SCF complexes, is required for the cellular transition from quiescence to the proliferative state. Mutations in this gene are associated with multiple endocrine neoplasia type IV (MEN4).

For Research Use Only

A Reliable Research Partner in Life Science and Medicine

Tel: 400-999-2100

Email: techsupport@elabscience.cn

Web: www.elabscience.cn