

SSSCA1 Polyclonal Antibody

Catalog Number:E-AB-52845

Note: Centrifuge before opening to ensure complete recovery of vial contents.


Description

Reactivity	Human, Mouse
Immunogen	Fusion protein of human SSSCA1
Host	Rabbit
Isotype	IgG
Purification	Affinity purification
Conjugation	Unconjugated
Formulation	PBS with 0.05% NaN ₃ and 40% Glycerol,pH7.4

Applications Recommended Dilution


WB	1:500-1:2000
IHC	1:50-1:300
ELISA	1:5000-1:10000

Data

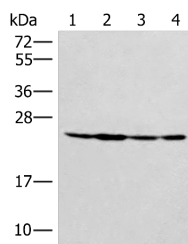


Western blot analysis of HeLa 231 and A549 cell
Mouse liver tissue lysates using SSSCA1 Polyclonal
Antibody at dilution of 1:400

Observed Mw:Refer to figures
Calculated Mw:21 kDa

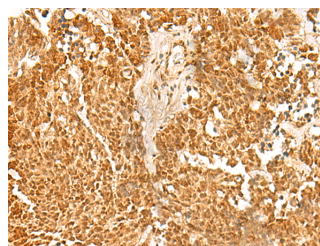

Immunohistochemistry of paraffin-embedded
Human lung cancer tissue using SSSCA1 Polyclonal
Antibody at dilution of 1:50(×200)

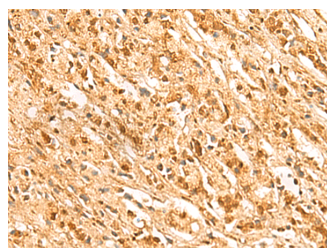

Immunohistochemistry of paraffin-embedded
Human prostate cancer tissue using SSSCA1
Polyclonal Antibody at dilution of 1:50(×200)

Preparation & Storage

Storage Store at -20°C. Avoid freeze / thaw cycles.

Background

For Research Use Only

A Reliable Research Partner in Life Science and Medicine

Tel: 400-999-2100

Email: techsupport@elabscience.cn

Web: www.elabscience.cn

SSSCA1 Polyclonal Antibody

Catalog Number:E-AB-52845



This antigen is recognized by a subset of anti-centromere antibodies from patients with scleroderma and/or Sjogren's syndrome. Subcellular localization has not yet been established.

For Research Use Only

A Reliable Research Partner in Life Science and Medicine

Tel: 400-999-2100

Email: techsupport@elabscience.cn

Web: www.elabscience.cn