

SARS-COV-2 NP Polyclonal Antibody(2019-nCoV)

Catalog Number:E-AB-V1015



Note: Centrifuge before opening to ensure complete recovery of vial contents.


Description

Reactivity	SARS-COV2
Immunogen	Recombinant 2019-nCoV Nucleocapsid Protein (His Tag)
Host	Rabbit
Isotype	IgG
Conjugation	Unconjugated
Formulation	0.2 µm filtered solution in PBS

Applications Recommended Dilution

WB	1:1000-1:5000_x000D_
ELISA:1:5000-1:1000	
0	

Data

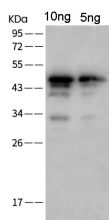

Western Blot analysis of SARS-CoV2-Spike protein using SARS-COV/SARS-COV-2 NP Polyclonal Antibody(2019-nCoV) at dilution of 1:2000

Preparation & Storage

Storage Store at -20°C. Avoid freeze / thaw cycles.

Background

Coronavirus N protein is required for coronavirus RNA synthesis, and has RNA chaperone activity that may be involved in template switch. Nucleocapsid protein is a most abundant protein of coronavirus. N protein packages the positive strand viral genome RNA into a helical ribonucleocapsid (RNP) and plays a fundamental role during virion assembly through its interactions with the viral genome and membrane protein M. Plays an important role in enhancing the efficiency of subgenomic viral RNA transcription as well as viral replication. Because of the conservation of N protein sequence and its strong immunogenicity, the N protein of coronavirus is chosen as a diagnostic tool.

For Research Use Only

A Reliable Research Partner in Life Science and Medicine

Tel: 400-999-2100

Email: techsupport@elabscience.cn

Web: www.elabscience.cn