

ERK 1/2 Polyclonal Antibody

Catalog Number:E-AB-31374

3 Publications

Note: Centrifuge before opening to ensure complete recovery of vial contents.


Description

Reactivity	Human,Mouse,Rat
Immunogen	Synthesized peptide derived from the C-terminal region of human ERK 1/2
Host	Rabbit
Isotype	IgG
Purification	Affinity purification
Conjugation	Unconjugated
Formulation	PBS with 0.02% sodium azide, 0.5% BSA and 50% glycerol, pH7.4

Applications Recommended Dilution

WB	1:500-1:2000
IHC	1:100-1:300
IF	1:50-1:200
ELISA	1:10000


Data

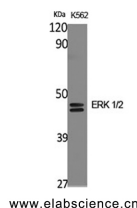


Western Blot analysis of K562 cells using ERK 1/2 Polyclonal Antibody at dilution of 1:2000.

Observed Mw:42,44kDa

Calculated Mw:43kDa

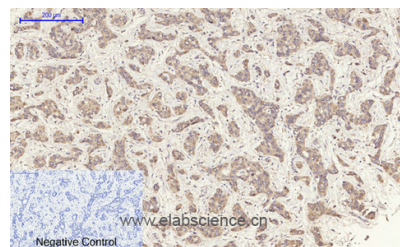

Immunohistochemistry of paraffin-embedded Human liver cancer tissue using ERK 1/2 Polyclonal Antibody at dilution of 1:200.

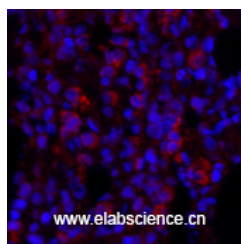

Immunofluorescence analysis of Rat lung tissue using ERK 1/2 Polyclonal Antibody at dilution of 1:200.

Preparation & Storage

Storage Store at -20°C. Avoid freeze / thaw cycles.

Background

Involved in both the initiation and regulation of meiosis, mitosis, and postmitotic functions in differentiated cells by

For Research Use Only

A Reliable Research Partner in Life Science and Medicine

Tel: 400-999-2100

Email: techsupport@elabscience.cn

Web: www.elabscience.cn

ERK 1/2 Polyclonal Antibody

Catalog Number: E-AB-31374

3 Publications



phosphorylating a number of transcription factors such as ELK1. Phosphorylates EIF4EBP1; required for initiation of translation. Phosphorylates microtubule-associated protein 2 (MAP2). Phosphorylates SPZ1 (By similarity).

Phosphorylates heat shock factor protein 4 (HSF4) and ARHGEF2. Acts as a transcriptional repressor. Binds to a [GC]AAA[GC] consensus sequence. Repress the expression of interferon gamma-induced genes. Seems to bind to the promoter of CCL5, DMP1, IFIH1, IFITM1, IRF7, IRF9, LAMP3, OAS1, OAS2, OAS3 and STAT1. Transcriptional activity is independent of kinase activity.

For Research Use Only

A Reliable Research Partner in Life Science and Medicine

Tel: 400-999-2100

Email: techsupport@elabscience.cn

Web: www.elabscience.cn