

BCL2 Polyclonal Antibody

Catalog Number: E-AB-15522

7 Publications



Note: Centrifuge before opening to ensure complete recovery of vial contents.

Description


Reactivity	Human
Immunogen	Synthetic peptide of human BCL2
Host	Rabbit
Isotype	IgG
Purification	Affinity purification
Conjugation	Unconjugated
Formulation	PBS with 0.05% sodium azide, 50% glycerol, PH7.3

Applications


Recommended Dilution

WB	1:500-1:2000
-----------	--------------


Data

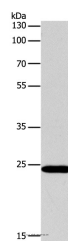

Western Blot analysis of Hela cell using BCL2 Polyclonal Antibody at dilution of 1:650
Calculated Mw:22kDa

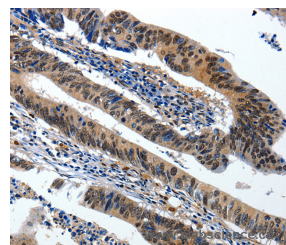

Immunohistochemistry of paraffin-embedded Human colorectal cancer using BCL2 Polyclonal Antibody at dilution of 1:20

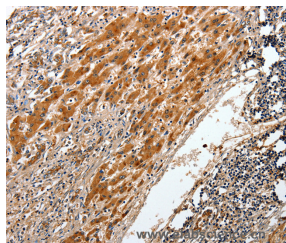

Immunohistochemistry of paraffin-embedded Human liver cancer using BCL2 Polyclonal Antibody at dilution of 1:20

Preparation & Storage

Storage Store at -20°C. Avoid freeze / thaw cycles.

Background

This gene encodes an integral outer mitochondrial membrane protein that blocks the apoptotic death of some cells such as lymphocytes. Constitutive expression of BCL2, such as in the case of translocation of BCL2 to Ig heavy chain locus, is thought to be the cause of follicular lymphoma. Two transcript variants, produced by alternate splicing, differ in their C-terminal ends.

For Research Use Only

A Reliable Research Partner in Life Science and Medicine

Tel: 400-999-2100

Email: techsupport@elabscience.cn

Web: www.elabscience.cn