

FOXO1 Polyclonal Antibody

Catalog Number:E-AB-70143



Note: Centrifuge before opening to ensure complete recovery of vial contents.


Description

Reactivity	Mouse,Rat
Immunogen	Recombinant protein corresponding to Mouse FOXO1
Host	Rabbit
Isotype	IgG
Purification	Affinity purification
Conjugation	Unconjugated
Formulation	PBS with 0.02% sodium azide,100 µg/ml BSA and 50% glycerol.

Applications Recommended Dilution

WB	1:500-1:2000
-----------	--------------

Data

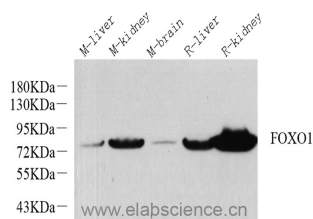

Western Blot analysis of various samples using FOXO1 Polyclonal Antibody at dilution of 1:1000.

Observed Mw:78kDa
Calculated Mw:78-82kDa

Preparation & Storage

Storage Store at -20°C. Avoid freeze / thaw cycles.

Background

FOXO1, also named as FOXO1A, FKHR and FKH1, is a member of the FOXO subfamily of Forkhead transcription factors. FOXO1 is a transcription factor which acts as a regulator of cell responses to oxidative stress. FOXO1 interacts with LRPPRC and SIRT1. In the presence of KIRT1, FOXO1 mediates down-regulation of cyclin D1 and up-regulation of CDKN1B levels which are required for cell transition from proliferative growth to quiescence. FOXO1 contains three predicted protein kinase B phosphorylation sites (Thr-24, Ser-256, and Ser-319) that are conserved in other FOXO proteins. The t(2;13) and the variant t(1;13) translocations generate PAX3/FKHR and PAX7/FKHR fusion proteins respectively. The resulting protein is a transcriptional activator. Defects in FOXO1 are a cause of rhabdomyosarcoma type 2 (RMS2).

For Research Use Only

A Reliable Research Partner in Life Science and Medicine

Tel: 400-999-2100

Email: techsupport@elabscience.cn

Web: www.elabscience.cn