

Caspase-3 Polyclonal Antibody

Catalog Number:E-AB-60017

8 Publications



Note: Centrifuge before opening to ensure complete recovery of vial contents.


Description

Reactivity	Human,Mouse,Rat
Immunogen	A synthetic peptide of human Caspase-3 (NP_004337.2).
Host	Rabbit
Isotype	IgG
Purification	Affinity purification
Conjugation	Unconjugated
Formulation	PBS with 0.02% sodium azide, 50% glycerol, pH7.3

Applications Recommended Dilution


WB	1:500-1:2000
IF	1:50-1:200

Data

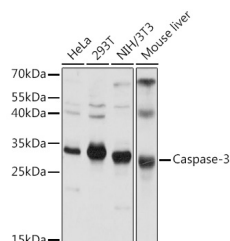

Western blot analysis of extracts of various cell lines using Caspase-3 Polyclonal Antibody at dilution of 1:1000.

Observed Mw:15 kDa/35 kDa
Calculated Mw:31 kDa

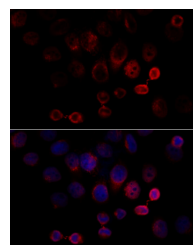

Immunofluorescence analysis of HeLa cells using Caspase-3 Polyclonal Antibody at dilution of 1:100 (40x lens). Blue: DAPI for nuclear staining.

Preparation & Storage

Storage Store at -20°C. Avoid freeze / thaw cycles.

Background

This gene encodes a protein which is a member of the cysteine-aspartic acid protease (caspase) family. Sequential activation of caspases plays a central role in the execution-phase of cell apoptosis. Caspases exist as inactive proenzymes which undergo proteolytic processing at conserved aspartic residues to produce two subunits, large and small, that dimerize to form the active enzyme. This protein cleaves and activates caspases 6, 7 and 9, and the protein itself is processed by caspases 8, 9 and 10. It is the predominant caspase involved in the cleavage of amyloid-beta 4A precursor protein, which is associated with neuronal death in Alzheimer's disease. Alternative splicing of this gene results in two transcript variants that encode the same protein.

For Research Use Only

A Reliable Research Partner in Life Science and Medicine

Tel: 400-999-2100

Email: techsupport@elabscience.cn

Web: www.elabscience.cn